

AUTOMATIC GRINDING & POLISHING MACHINE



DESCRIPTION

Model HUTOPOL-2000XP is a kind of double-disc pneumatic type automatic metallographic specimen grinding and polishing machine equipped with an automatic polishing head and step-less speed changing. It features advanced performance and high degree automation and can finish the whole process of rough grinding, fine grinding, rough polishing and finishing polishing automatically. It is the necessary metallographic specimen-preparing instrument for using in factories, scientific research institutes and laboratories of colleges.

Registered Office: UGF-12 Ansal Plaza-1 Chiranjeev Vihar Ghaziabad-201002 (U.P.) INDIA Workshop: Khasra No - 63 Village Bayana, Near Wave City, NH 24, Ghaziabad -201015 (U.P.) INDIA

Phone: +91-8750305291

Email : info@hsmleindia.in / sales@hsmleindia.com
Website : www.hopesmluckyenterprises.in / www.hsmleindia.com (Under Construction)



ISO 9001: 2015

FEATURES:

- Individual samples are loaded independently.
- Central force loading of the sample disc, fix six samples at one time, complete all grinding
 and polishing processes, and ensure the complete plane of each sample by the grinding.
- Grinding & polishing disc and sample disc: set and run all grinding and polishing parameters, such as rotating speed, grinding and polishing time, rotation direction, water valve on/off, etc., which are automatically saved and easy to call.
- Touch screen interface: it is convenient to set grinding and polishing parameters, intuitive status display and simple operation.
- The rotation direction of both grinding & polishing disc and sample disc can be switched, and the rotation speed changing is stepless.

 Working modes:
- (1) Fully automatic mode:

According to the sample materials or the habits of users, 30 sets of processes (processes) can be set and called, and each process can contain 10 process parameters (for each process, for a certain grinding or polishing process, process parameters: grinding disc and grinding head rotation speed, grinding and polishing time...)

(2) Single operation mode:

According to each grinding and polishing process, 30 process parameters can be set and called:

grinding disc and grinding head rotation speed, grinding, and polishing time, water on-

(3) Manual mode:

Operate a function of the equipment separately

- Four channel automatic dispenser (optional):
- (1) The main machine of the grinding and polishing machine controls the four channel dispenser. The type and speed of dispenser are controlled by the grinding and polishing machine.
- (2) There are three working modes of the dispenser: manual, online and automatic.
- (3) Connect the dispenser with the grinding & polishing machine and drip liquid according to the parameters set by the grinding and polishing machine. Each mode can set dripping speed and dripping time independently.

Registered Office: UGF-12 Ansal Plaza-1 Chiranjeev Vihar Ghaziabad-201002 (U.P.) INDIA Workshop: Khasra No - 63 Village Bayana, Near Wave City, NH 24, Ghaziabad -201015 (U.P.) INDIA


Phone: +91-8750305291

Email: info@hsmleindia.in / sales@hsmleindia.com

Website: www.hopesmluckyenterprises.in / www.hsmleindia.com (Under Construction)



- The USB flash disk stores 300 processes, which can be imported and exported.
- Electromagnetic automatic locking of grinding head, replacing manual grinding head locking wrench, which is more convenient.
- By easily replacing the magnetic anti-sticky disc, all processes such as rough and fine grinding, rough and fine polishing of various samples can be completed. One disc is equivalent to several discs.
- Mobile terminal remote control operation (optional).

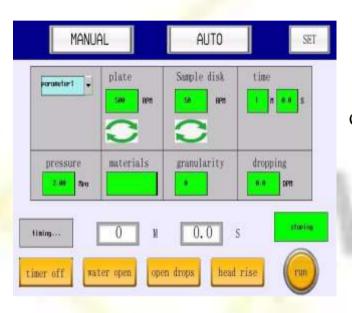

Operation interface





Electromagnetic automatic lock & remote

Registered Office: UGF-12 Ansal Plaza-1 Chiranjeev Vihar Ghaziabad-201002 (U.P.) INDIA Workshop: Khasra No - 63 Village Bayana, Near Wave City, NH 24, Ghaziabad -201015 (U.P.) INDIA

Phone: +91- 8750305291

Email: info@hsmleindia.in / sales@hsmleindia.com

Website: www.hopesmluckyenterprises.in / www.hsmleindia.com (Under Construction)

ISO 9001: 2015

TECHNICAL SPECIFICATIONS:

Model No.	HUTOPOL-2000XP
Grinding & polishing disc	254mm (optional 230mm, 203mm) with magnetic disc
diameter	
Grinding disc rotation speed	50~1000 rpm (Forward and reverse rotation direction can
	be switched)
Sample disc speed	30~200 rpm (Forward and reverse rotation direction can
	be switched)
Manual mode	30 sets of parameters could be selected, each one can be
	set and called separately
Automatic mode	30 sets of processes, each set of processes contains 10-step
	parameters (from step 1 to step 10 automatically, each
	step has an independent set of parameters)
Sample holder	Ø30mm, 6-hole (optional 20mm, 6-hole, 40mm/50mm 3-
	hole, or others can be customized)
Grinding head-locking mode	Electromagnetic automatic lock
Working air pressure	Not less than 0.6 MPa
Motor power	0.75Kw
Power supply	AC 220V, 50/60Hz, Single phase, 2KW
Dimensions & Weight	760 x 730 x 700mm & 74kg

STANDARD ACCESSORIES:

Item Name	Q'ty	Item Name	Q'ty
HUTOPOL 2000XP	1 set	Diamond suspension polishing liquid, 6	1 Pc
		microns, 200ml	
One side sticky Abrasive paper	2 Pc	One side sticky polishing cloth (250mm),	1 Pc
(250mm) 180#		flocking (brown)	
One side sticky Abrasive paper	2 Pc	One side sticky polishing cloth (250mm),	1 Pc
(250mm) 800#		canvas	
One side sticky Abrasive paper	2 Pc	Anti-sticky disc (250mm)	2 Pc
(250mm) 1000#			
Waterproof ring	2 Pc	Magnetic working plate (250mm)	2 Pc
Input-water pipe 8mm	1 Pc	Power cable (Aviation plug joint)	1Pc
Output-water pipe 25mm	1 Pc	Data cable (Aviation plug joint)	1 Pc
Hose hoops 18-32	1 Pc	Triangle valve Φ12.7mm	1 Pc
Hook	1 Pc	Teflon tape	1 Pc
Rubber pad with strainer Φ12.7mm	1 Pc	Fast twist joint SNZ08-04	1 Pc
Touch pen	1 Pc	Filter piece Φ12.7mm	4 Pc
Technical documents: Machine manual, certificate, packing list 1 share for each			

Registered Office: UGF-12 Ansal Plaza-1 Chiranjeev Vihar Ghaziabad-201002 (U.P.) INDIA Workshop: Khasra No - 63 Village Bayana, Near Wave City, NH 24, Ghaziabad -201015 (U.P.) INDIA

Phone: +91- 8750305291

Email: info@hsmleindia.in / sales@hsmleindia.com
Website: www.hopesmluckyenterprises.in / www.hsmleindia.com (Under Construction)



OPTIONAL ACCESSORIES:

Air pump (optional)				
Air pipe 6mm x 3m (optional)				
Huawei M6 tablet PC (10.2", 4GB+64GB Wi-Fi)				
A 4-channel polishing liquid dispen <mark>ser</mark>				
Metallographic	GSC250A	• 100 piec <mark>es /p</mark> ackage		
grinding paper		Material: black carborundum		
		• Sand size: P80# 120# 180# 320# 600# 800# 1000#		
		1200# 1500# 2000#		
		• Size: Φ250mm		
		Waterproof and one side sticky, can be attached		
		on the grinding plate directly		
Anti-sticky	MAGNOMET	Anti-sticky disc: 5 pieces/box		
disc		• Diameter: Φ <mark>250mm</mark>		
		 Surface coating with Special anti-sticky layer, 		
. 30	31 //	can attach th <mark>e stick</mark> y paper and polishing cloth		
- 10	off No.	directly, and easily replace without sticky mark.		
	11	Can be reused		
Polishing cloth	FT250A	Polishing cloth 10 pieces/package		
		Velvet fabric(long floss) / canvas (short floss)/		
		Woolen cloth (short floss) / silk cloth (no floss)		
		• Φ250mm		
		 One side sticky, can be attached directly on the 		
		grinding plate		
Diamond	DSU	 Diamond grinding and polishing liquid, 		
grinding and		500ml/bottle		
polishing		Diamond suspending in the liquid, can be		
liquid		dropped automatically or watering pot,		
		Uniformly distribution		
		• Size: 0.25, 0.5, 1, 1.5, 2.5, 3, 3.5, 5, 6, 7, 9, 10, 14,		
		15, 20, 28, 40micron, for different material		

Registered Office: UGF-12 Ansal Plaza-1 Chiranjeev Vihar Ghaziabad-201002 (U.P.) INDIA Workshop: Khasra No - 63 Village Bayana, Near Wave City, NH 24, Ghaziabad -201015 (U.P.) INDIA

Phone: +91-8750305291

Email: info@hsmleindia.in / sales@hsmleindia.com
Website: www.hopesmluckyenterprises.in / www.hsmleindia.com (Under Construction)