

AUTOMATIC METALLOGRAPHIC SPECIMEN GRINDING & POLISHING MACHINE



FEATURES:

 Multi-usage, one machine to complete the whole process of metallographic rough grinding, fine grinding, rough polishing and fine polishing.

ISO 9001: 2015

- Single-point pneumatic holding, no need for sample clamping.
- Six pieces of Φ30mm samples can be polished simultaneously.
- PLC independent control for the grinding disc and grinding head. The grinding & polishing parameters like rotation speed, grinding & polishing time, rotation direction, water valve on/off etc. can be set and automatically saved, easy to call.
- Large touchscreen interface, convenient parameter setting, intuitive state display and easy operation.
- Stepless speed changing for both grinding disc and grinding head. The rotation direction is switchable between FWD & REV.
- PLC control for water supply and grinding material dispenser.

 Working mode:
- (1) Automatic mode: According to the sample materials user's habits, 30 groups of processes can be set and called. Each process can contain 10 steps of process parameters (each process is corresponding to certain grinding or polishing process. Process parameters: rotation speed of grinding disc and polishing head, grinding and polishing time, etc.)
- (2) Manual mode: According to each grinding and polishing process, 30 groups of process parameters can be set and called: grinding disc and polishing head rotation speed, grinding and polishing time, water supply on/off etc.
- (3) Setting mode: Set certain function of equipment separately.

Registered Office: UGF-12 Ansal Plaza-1 Chiranjeev Vihar Ghaziabad-201002 (U.P.) INDIA Workshop: Khasra No - 63 Village Bayana, Near Wave City, NH 24, Ghaziabad -201015 (U.P.) INDIA

Phone: +91-8750305291

Email : info@hsmleindia.in / sales@hsmleindia.com
Website : www.hopesmluckyenterprises.in / www.hsmleindia.com (Under Construction)




HSMLE India Corporation

- ISO 9001: 2015
- Magnetic anti-sticky disc can be replaced easily and it can finish the rough grinding, fine grinding, rough polishing and fine polishing for various samples continuously.
- The water supply can be controlled by electromagnetic valve. There is automatic water cleaning device under grinding disc and it can clean the inside wall automatically.
- Grinding head automatic electromagnetic lock, no need for manual locking handle.
- Multi-language interface is available.
- High accuracy, stable running, low noise.
- Grinding & polishing dispenser is optional.
- Operation by mobile remote control device is optional.

EATURE:





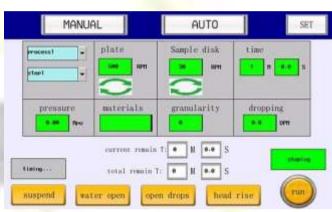

ECHNICAL SPECIFICATIONS

Model No.	HUTOPOL-1000
Magnetic grinding disc	250mm (optional 200mm, 230mm)
diameter	
Sample holder	6-hole, Φ30mm (Φ20mm, or 3-hole, Φ40-Φ50mm is
	optional)
Grinding disc rotation speed	100-1000 rpm, FWD & REV rotation directions can be
	switched
Grinding head rotation speed	30-200 rpm, FWD & REV rotation directions can be
	switched
Manual mode	Manual mode
Automatic mode	30 groups of process parameters, and each group of
	process includes 10 steps of procedure parameters
	(automatic processing from step 1 to step 10, each step
	has its independent group of procedure parameters).

Registered Office: UGF-12 Ansal Plaza-1 Chiranjeev Vihar Ghaziabad-201002 (U.P.) INDIA Workshop: Khasra No - 63 Village Bayana, Near Wave City, NH 24, Ghaziabad -201015 (U.P.) INDIA

Phone: +91-8750305291

Email: info@hsmleindia.in / sales@hsmleindia.com

Website: www.hopesmluckyenterprises.in / www.hsmleindia.com (Under Construction)



HSMLE India Corporation

IS	n .	a	n	n	1	2	n	4	5
10		21	2.4	•			v	u	•

Grinding head-locking	Electromagnetic automatic locking	
method		
Air source pressure	no less than 0.6MPa	
Power supply	220V, 50/60Hz, single phase, 2KW	
Motor power	750W	
Overall Dimensions	760x 470 x 700 mm	
Net weight	80Kg	

STANDARD ACCESSORIES:

Item Name	Q'ty	Item Name	Q'ty	
HUTOPOL-1000	1 set	Diamond suspension polishing	1	
		liquid, 6 micron, 200ml	bottle	
One side sticky Abrasive paper	2 pc	One side sticky polishing cloth	1 pc	
(250mm) 180#		(250mm), flocking (brown)		
One side sticky Abrasive paper	2 pc	One side sticky polishing cloth	1 pc	
(250mm) 800#		(25 <mark>0mm),</mark> canvas		
One side sticky Abrasive paper	2 pc	Anti-sticky disc (250mm)	1 pc	
(250mm) 1000#				
Waterproof ring	1 pc	Magnetic working plate (250mm)	1 pc	
Input-water pipe 8mm	1 pc	Power cable	1 pc	
		(Aviation plug joint)		
Output-water pipe 25mm	1 pc	Data cable (Aviation plug joint)	1 pc	
Hose hoops 18-32	1 pc	Triangle valve Φ12.7mm	1 pc	
Hook	1 pc	Teflon tape	1 pc	
Rubber pad with strainer Φ12.7mm	1 pc	Fast twist joint SNZ08-04	1 pc	
Touch pen	1 pc	Filter piece Φ12.7mm	4 pcs	
Technical documents: Machine manual, certificate, packing list 1 share for each				



Registered Office: UGF-12 Ansal Plaza-1 Chiranjeev Vihar Ghaziabad-201002 (U.P.) INDIA Workshop: Khasra No - 63 Village Bayana, Near Wave City, NH 24, Ghaziabad -201015 (U.P.) INDIA

Phone: +91-8750305291

Email: info@hsmleindia.in / sales@hsmleindia.com
Website: www.hopesmluckyenterprises.in / www.hsmleindia.com (Under Construction)



OPTIONAL ACCESSORIES:

Air pump (optional)				
Air pipe 6mmx3m (optional)				
Tablet PC (10.2", 4GB+64	4GB WiFi)			
Metallographic grinding paper	GSC250A	 100 pieces /package Material: black carborundum Sand size: P80# 120# 180# 320# 600# 800# 1000# 1200# 1500# 2000# Size: Φ250mm Waterproof and one side sticky, can be attached on the grinding plate directly 		
Anti-sticky disc	MAGNOMET	 Anti-sticky disc: 5 pieces/box Diameter: Φ250mm Surface coating with Special anti-sticky layer, can attach the sticky paper and polishing cloth directly, and easily replace without sticky mark. Can be reused 		
Polishing cloth	FT250A	 Polishing cloth 10 pieces/package Velvet fabric(long floss) / canvas (short floss)/ Woolen cloth (short floss) / silk cloth (no floss) Ф250mm One side sticky, can be attached directly on the grinding plate 		
Diamond grinding and polishing liquid	DSU	 Diamond grinding and polishing liquid, 500ml/bottle Diamond suspending in the liquid, can be dropped automatically or watering pot, Uniformly distribution Size: 0.25, 0.5, 1, 1.5, 2.5, 3, 3.5, 5, 6, 7, 9, 10, 14, 15, 20, 28, 40micron, for different material 		

Registered Office: UGF-12 Ansal Plaza-1 Chiranjeev Vihar Ghaziabad-201002 (U.P.) INDIA Workshop: Khasra No - 63 Village Bayana, Near Wave City, NH 24, Ghaziabad -201015 (U.P.) INDIA

Phone: +91-8750305291

Email: info@hsmleindia.in / sales@hsmleindia.com
Website: www.hopesmluckyenterprises.in / www.hsmleindia.com (Under Construction)