

SINGLE-DISC AUTOMATIC GRINDING AND POLISHING MACHINE



DESCRIPTION:

ISO 9001 : 2015

Model HUTOPOL-1000XP is a kind of single disc pneumatic type automatic metallographic specimen grinding and polishing machine equipped with automatic polishing head and step-less speed changing. It features advanced performance and high degree automation and can finish the whole process of rough grinding, fine grinding, rough polishing and finishing polishing automatically. It is the metallographic necessary specimen-preparing instrument for using in factories, scientific research institutes and laboratories of colleges.

FEATURES:

- Individual samples are loaded independently.
- Central force loading of the sample disc, fix six samples at one time, complete all
 grinding and polishing processes, and ensure the complete plane of each sample by
 the grinding.
- Grinding & polishing disc and sample disc: set and run all grinding and polishing parameters, such as rotating speed, grinding and polishing time, rotation direction, water valve on/off, etc., which are automatically saved and easy to call.
- Touch screen interface: it is convenient to set grinding and polishing parameters, intuitive status display and simple operation.

Registered Office: UGF-12 Ansal Plaza-1 Chiranjeev Vihar Ghaziabad-201002 (U.P.) INDIA Workshop: Khasra No - 63 Village Bayana, Near Wave City, NH 24, Ghaziabad -201015 (U.P.) INDIA

Phone: +91-8750305291

Email: info@hsmleindia.in / sales@hsmleindia.com

Website: www.hopesmluckyenterprises.in / www.hsmleindia.com (Under Construction)

ISO 9001: 2015

• The rotation direction of both grinding & polishing disc and sample disc can be switched, and the rotation speed changing is stepless.

Working modes:

(1) Fully automatic mode:

According to the sample materials or the habits of users, 30 sets of processes (processes) can be set and called, and each process can contain 10 process parameters (for each process, for a certain grinding or polishing process, process parameters: grinding disc and grinding head rotation speed, grinding and polishing time...)

(2) Single operation mode:

According to each grinding and polishing process, 30 process parameters can be set and called: grinding disc and grinding head rotation speed, grinding and polishing time, water on-off

(3) Manual mode:

Operate a function of the equipment separately


- Four-channel automatic dispenser (optional):
- (1) The main machine of the grinding and polishing machine controls the four-channel dispenser. The type and speed of the dispenser are controlled by the grinding and polishing machine.
- (2) There are three working modes of the dispenser: manual, online and automatic.
- (3) Connect the dispenser with the grinding & polishing machine and drip liquid according to the parameters set by the grinding and polishing machine. Each mode can set dripping speed and dripping time independently.
- The USB flash disk stores 300 processes, which can be imported and exported.
- Electromagnetic automatic locking of grinding head, replacing manual grinding head locking wrench, which is more convenient.
- By easily replacing the magnetic anti sticky disc, all processes such as rough and fine grinding, rough and fine polishing of various samples can be completed. One disc is equivalent to several discs.
- Mobile terminal remote control operation (optional).

Registered Office: UGF-12 Ansal Plaza-1 Chiranjeev Vihar Ghaziabad-201002 (U.P.) INDIA Workshop: Khasra No - 63 Village Bayana, Near Wave City, NH 24, Ghaziabad -201015 (U.P.) INDIA

Phone: +91- 8750305291

Email: info@hsmleindia.in / sales@hsmleindia.com

Website: www.hopesmluckyenterprises.in / www.hsmleindia.com (Under Construction)

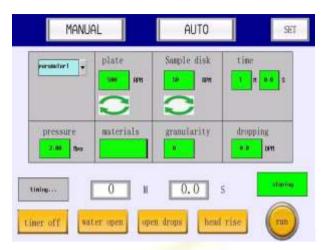





Operation interface

Electromagnetic automatic lock & remote

TECHNICAL SPECIFICATIONS:

Model No.	HUTOPOL-1000XP		
Grinding & polishing	254mm (optional 230mm, 203mm) with magnetic disc		
disc diameter	AND DOCUMENTS AND DESCRIPTION OF THE PARTY O		
Grinding disc rotation	50~1000 rpm (Forward and reverse rotation direction can be		
speed	switched)		
Sample disc speed	30~200 rpm (Forward and reverse rotation direction can be		
	switched)		
Manual mode	30 sets of parameters could be selected, each one can be set		
	and called separately		
Automatic mode	30 sets of processes, each set of process contains 10-step		
	parameters (from step 1 to step 10 automatically, each step		
	has an independent setting of parameters)		
Sample holder	Ø30mm, 6-hole (optional 20mm, 6-hole, 40mm/50mm 3-hole,		
	or others can be customized)		
Grinding head-	ad- Electromagnetic automatic lock		
locking mode			
Working air pressure	Not less than 0.6 MPa		
Motor power	0.75Kw		
Power supply	AC 220V, 50/60Hz, Single phase, 2KW		
Dimensions & Weight	760 x 730 x 700mm & 63kg		

Registered Office: UGF-12 Ansal Plaza-1 Chiranjeev Vihar Ghaziabad-201002 (U.P.) INDIA Workshop: Khasra No - 63 Village Bayana, Near Wave City, NH 24, Ghaziabad -201015 (U.P.) INDIA

Phone: +91-8750305291

Email: info@hsmleindia.in / sales@hsmleindia.com
Website: www.hopesmluckyenterprises.in / www.hsmleindia.com (Under Construction)





STANDARD ACCESSORIES:

Item Name	Q'ty	Item Name	Q'ty
HUTOPOL 1000XP	1	Diamond suspension polishing	1 Pc
	set	liquid, 6 micron 200ml	
One side sticky Abrasive paper	2 Pc	One side sticky polishing cloth	1 Pc
(250mm) 180#		(250mm), flocking (brown)	
One side sticky Abrasive paper	2 Pc	One side sticky polishing cloth	1 Pc
(250mm) 800#		(250 <mark>mm), c</mark> anvas	
One side sticky Abras <mark>ive pa</mark> per	2 Pc	Anti-sticky disc (250mm)	2 Pc
(250mm) 1000#			
Waterproof ring	2 Pc	Magnetic working plate (250mm)	2 Pc
Input-wate <mark>r pipe</mark> 8mm	1 Pc	Power cable (Aviation plug joint)	1Pc
Outp <mark>ut-wate</mark> r pipe 25mm	1 Pc	Data cable (Aviation plug joint)	1 Pc
Hose hoops 18-32	1 Pc	Triangle valve Φ12.7mm	1 Pc
Hook	1 Pc	Tef <mark>lon tap</mark> e	1 Pc
Rubber pad with strainer	1 Pc	Fast twist joint SNZ08-04	1 Pc
Φ12.7mm	- 675		
Touch pen	1 Pc	Filter p <mark>iec</mark> e Φ12.7mm	4 Pc
Technical documents: Machine m	anual, c	certificate <mark>, pac</mark> king list 1 share for each	

ISO 9001: 2015

Registered Office: UGF-12 Ansal Plaza-1 Chiranjeev Vihar Ghaziabad-201002 (U.P.) INDIA Workshop: Khasra No - 63 Village Bayana, Near Wave City, NH 24, Ghaziabad -201015 (U.P.) INDIA

Phone: +91-8750305291

Email: info@hsmleindia.in / sales@hsmleindia.com
Website: www.hopesmluckyenterprises.in / www.hsmleindia.com (Under Construction)



OPTIONAL ACCESSORIES:

Air pump (optio	onal)				
Air pipe 6mm x	Air pipe 6mm x 3m (optional)				
Huawei M6 tablet PC (10.2", 4GB+64GB WiFi)					
4-channel polishing liquid dispenser					
Metallographic	GSC250A	• 100 pieces /package			
grinding paper		Material: black carborundum			
		• Sand size: P80# 120# 180# 320# 600# 800# 1000#			
		1200# 1500# 2000#			
		• Size: Φ250mm			
		Waterproof and one side sticky, can be attached on			
		the grinding plate directly			
Anti-sticky	MAGNOMET	Anti-sticky disc: 5 pieces/box			
disc		• Diameter: Φ250mm			
		 Surface coating with Special anti-sticky layer, can 			
		attach the stick <mark>y pape</mark> r and polishing cloth directly,			
		and easily replace without sticky mark.			
		Can be reused			
Polishing cloth	FT250A	 Polishing cloth 10 pieces/package 			
		 Velvet fabric(long floss) / canvas (short floss)/ 			
		Woolen cloth (short floss) / silk cloth (no floss)			
		• Φ250mm			
		 One side sticky, can be attached directly on the 			
		grinding plate			
Diamond	DSU	 Diamond grinding and polishing liquid, 			
grinding and		500ml/bottle			
polishing		 Diamond suspending in the liquid, can be dropped 			
liquid		automatically or watering pot, Uniformly			
		distribution			
		• Size: 0.25, 0.5, 1, 1.5, 2.5, 3, 3.5, 5, 6, 7, 9, 10, 14, 15,			
		20, 28, 40micron, for different material			

Registered Office: UGF-12 Ansal Plaza-1 Chiranjeev Vihar Ghaziabad-201002 (U.P.) INDIA Workshop: Khasra No - 63 Village Bayana, Near Wave City, NH 24, Ghaziabad -201015 (U.P.) INDIA

Phone: +91-8750305291

Email: info@hsmleindia.in / sales@hsmleindia.com
Website: www.hopesmluckyenterprises.in / www.hsmleindia.com (Under Construction)